



MRI Core DX Extended Top

51817



BUILT FOR MRI. ENGINEERED FOR EFFICIENCY.

MRI-conditional storage cart with non-metal drawers, extended top work surface, and side compartments for safe, efficient organization.

MR CONDITIONAL DESIGN

- Three drawers, one sliding door compartment, six side compartments (2 dispensing, 4 enclosed), extended top surface (28 $\frac{7}{8}$ " x 19 $\frac{5}{8}$ ").

EXTENDED WORK SURFACE

- Extended top adds more workspace without increasing footprint.

SAFE FUNCTIONAL MATERIALS

- High-density polyethylene (HDPE), aluminum composite drawers, stainless steel fasteners.

BUILT FOR MRI FUNCTION

- MRI-conditional storage and prep surface in one cart
- Efficient access and mobility in metal-sensitive areas.

LIFETIME GUARANTEE

- 100% Lifetime Warranty, Made in USA



ASK ABOUT OUR DEMO PROGRAM

[Try it out at your facility.](#)



Sliding Door

- Secure items in the bulk storage area with the sliding door.



Non-Metal Slides

- Extend 75% and close smoothly with a durable plastic detent.



Dispensing Shelf

- Easily access gloves, kimwipes or other supplies.



Easy-Grab Handles

- Ergonomic, lightweight, and easy to maneuver.



10991 N Airworld Drive, Kansas City, MO, USA | 1 (800) 874-7768 | Sales@TrippNT.com

MRI Core DX with Extended Top


51817




Premium high density polyethylene surfaces.



High quality metal-free drawer slide with detent!

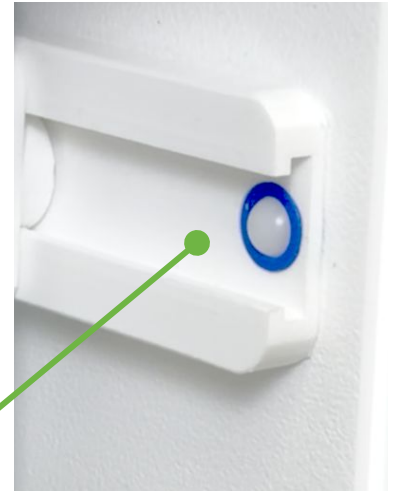

See Drawer In Action!

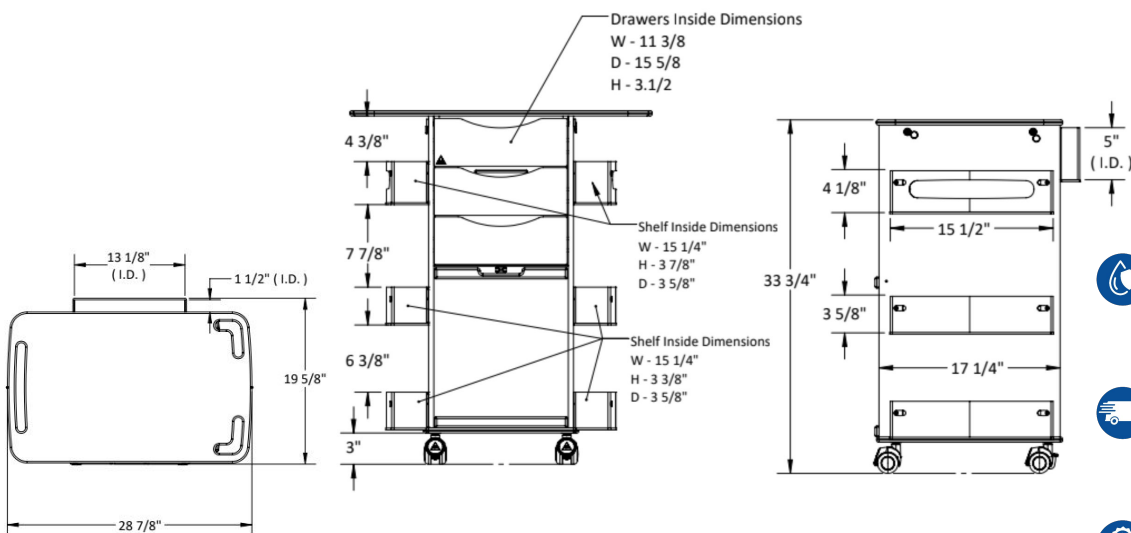

- Corrosion Resistant
- Short Lead Times
- Ships Assembled
- Made in the U.S.A
- 100% Lifetime Guarantee

ACCESSORIES

